

Capp *Rondo*


TUBE ROLLER

Capp *Rondo* Tube Roller ensures efficient mixing of blood samples, viscous substances and liquid-solid suspensions. The programmable pulse mode feature with Run-Pause-Run motion and variable time setting (30 to 99 seconds) give the user additional flexibility and improve the mixing process. In addition to the rolling motion, the rollers have a rocking movement which ensures optimal mixing throughout. Equipped with a brushless DC motor, **Capp *Rondo* Tube Roller** provides quiet and maintenance-free operation, highly suitable for extended runs.

Highlights:

Improved mixing performance:

- In addition to the rolling motion, the rollers of the **Capp *Rondo* Tube Roller** rock by lowering and raising, which improves the mixing process.
- Programmable momentary pulse mode ensures throughout optimal mixing.

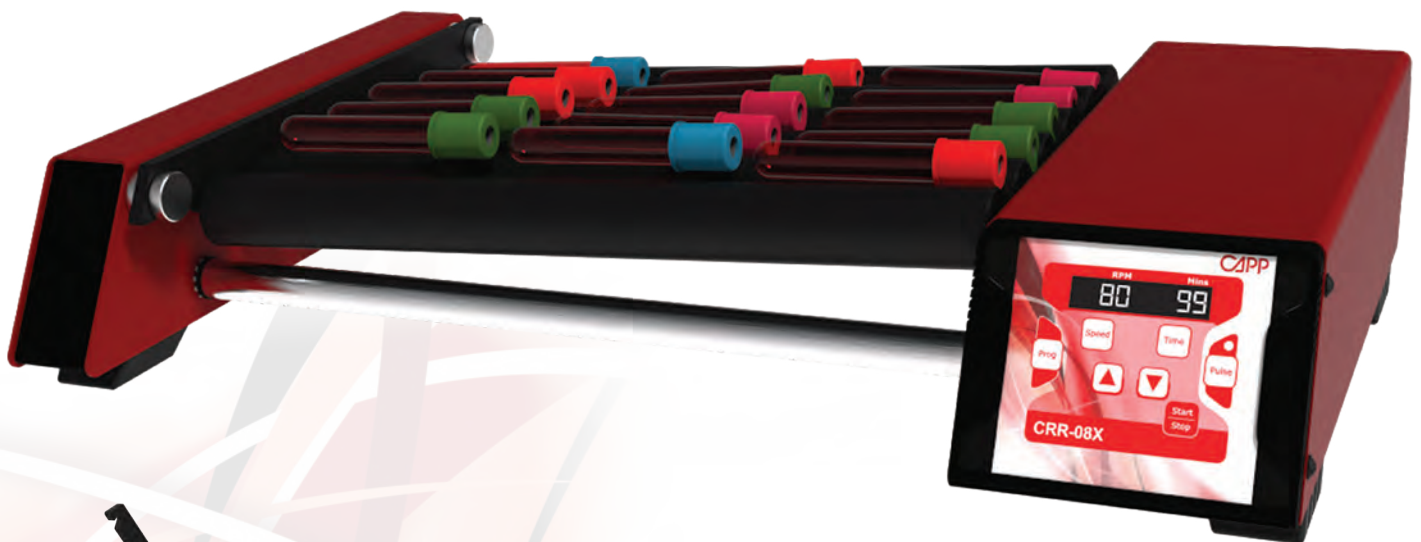

Removable Rollers



Stackable System



Capp Rondo Tube Roller can accommodate up to 15 blood tubes of 7mL. The user interface panel is highly intuitive and displays both the time and the speed.

Capp Rondo Tube Roller operates at the speed range from 10 up to 80 rpm. The user can program up to 9 customized protocols.

Focusing on increased user comfort and the low contamination risk, **Capp Rondo Tube Roller** operates with easily removable and autoclavable rollers. The stackable design and relatively small footprint will optimize your bench space. **Capp Rondo Tube Roller** is also suitable for use in incubators with temperatures up to 40°C and humidity up to 80%.



Highlights:

- Suitable for incubators and highly humid environments.
- Tool free detachable and autoclavable rollers for increased user comfort and lower contamination risk.
- By adjusting the rollers space, you can accommodate multiple sizes of tubes and bottles.
- Stackable design with a small footprint for optimal bench space use.

Specifications:	Capp Rondo Tube Roller
Motor type	Brushless DC motor
Motion	Rocking and rolling
Number of tubes	up to 15 x 7ml blood tubes
Speed range	10 to 80 rpm
Time range	1 to 99 mins
Programmable range	User can set and save up to 9 user-defined programs (protocols)
Operation type	Continuous and timed operation
Removable rollers	Yes, tool free removal and fitment
Pulse mode	Yes (Run-Pause-Run mode)
Tilt angle	+/-3°
Protection class	IP21
Weight	6kg
Dimensions (W x D x H)	502mm x 92.5mm x 301mm
Warranty	2 years

Ordering Info

Cat No.	Description
CRR-08X	CappRondo Tube Roller 80 rpm